




BEARING FRICTION APPARATUS

HFN6

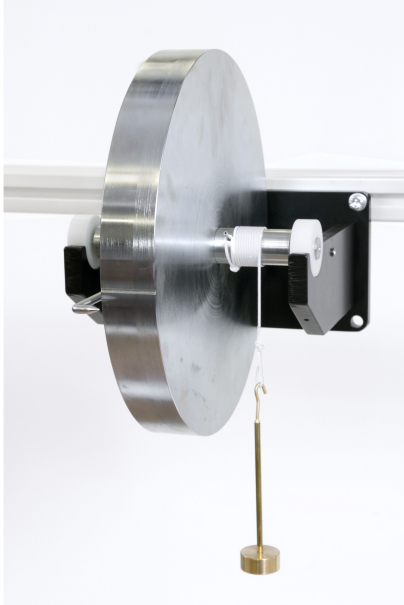

Year 2
study

Features

- Sturdy Wall mounted apparatus
- Can be mounted to HST1, HST100, HAC10, HAC14, HVT12f frames
- Compares Sliding and Roller Bearing efficiency
- Four solid bearings supplied as standard
- Can be used for Inertia and energy transfer experiments
- Customer specific bearings can be used
- Full set of weights and hanger supplied

Description

The apparatus consists of a flywheel on a horizontal shaft carried in a pair of similar bearings, the flywheel being used to even out small variations in friction. The shaft is of mild steel, and interchangeable bearings in a range of materials are provided. This experiment provides an opportunity to study journal bearing friction over a range of commonly used materials including a self lubricating sintered bronze. For comparison a ball or roller bearing is available to demonstrate the advantages of rolling rather than sliding contact. Torque is applied by a loaded cord wrapped round the flywheel shaft. A set of calibrated weights and Load hanger is supplied for loading the flywheel.

Related Laws/Applications

- Automotive
- Friction
- Machines
- Bearing
- Wear
- Torque
- Slip

Learning capabilities

- To compare the frictional losses of bearings by measuring the coefficient of sliding friction between pairs of materials
- To measure the reduced losses when a ball or roller bearing is used

Technical Specification

- Flywheel: Ø300 x 40mm thick; approximately 24kg
- Journal Bearing materials: Cast Iron, Tufnel, Sintered bronze, ball bearings

Recommended Ancillaries

- HAC10
- HAC14
- HST1
- HST100
- HVT12f

What's in the Box?

- 1 x HFN6 Flywheel
- 1 x Wall Bracket Assembly
- 2 x Load hanger
- 2 x Cord assemblies
- 1 x 5m Spare cord
- 2 x Tufnol bearing
- 2 x Sintered bronze bearing
- 2 x Cast iron bearing
- 2 x Ball bearing
- 10 x 0.1N; 10 x 0.2N; 1 x 1N; 2 x 2N; 1 x 5N; 2 x 10N
- Instruction manual
- Packing list
- Test sheet

You might also like

- HFN6A

Weights & Dimensions

- Weight: approximately 30kg
- Length: 200mm
- Width: 300mm
- Height: 400mm

Essential Services

- Rigid Vertical Support

Ordering information

To order this product, please call PA Hilton quoting the following code:
HFN6

All brand and/or product names are trademarks of their respective owners. Specifications and external appearance are subject to change without notice. The colour of the actual product may vary from the colour shown in the brochure.

Copyright © 2018 P.A. Hilton Limited. All rights reserved. This technical leaflet, its contents and/or layout may not be modified and/or adapted, copied in part or in whole and/or incorporated into other works without the prior written permission of P. A. Hilton Limited. Hi-Tech Education is a registered trade mark of P. A. Hilton Limited.

COUNTRY OF ORIGIN - UK WARRANTY PERIOD - 5 YEARS